

Emperor

EmFilm S Shaped Film Dressing

Made from latex-free, semi-permeable polyurethane film membrane and coated with a hypoallergenic acrylic adhesive, EmFilm S Framed Dressings are designed to be breathable with good oxygen and moisture vapour exchange.

Waterproof and impermeable to liquids, bacteria and viruses, the risk of infection is reduced while also enabling effective wound monitoring through a transparent film.



Protection of newly forming tissue



Fixation of primary dressings



Provision of a 'moist wound healing environment'



The patient can shower with the dressing in place



Prevention of contamination by dirt or bacteria, thus a reduced risk of wound infection



Rounded edges and consistent hypoallergenic adhesive help prevent dressing damage and early removal, reducing the number of dressing changes needed.

Call: 0800 008 6646

E-mail: woundcare@rociallehealthcare.com

Visit: www.emperorwound.com

 Rociale
Healthcare

Emperor is a product brand of Rociale Healthcare

Size	Product Code	NHSSC Code	PIP Code	Drug Tariff Price
5cm x 7cm	537708ZD	ELW85031	424-3705	£0.45
10cm x 12cm	537709ZD	ELW85044	424-3697	£1.04
10cm x 25cm	537710ZD	ELW85043	424-3671	£2.10
15cm x 20cm	537711ZD	ELW85041	424-3689	£2.04

For **FREE** product samples and more information, please get in touch


EmFilm S Application Guide

Preparation:

1. Ensure that the skin surrounding the wound is clean, dry, and free from grease, soaps or detergents
2. Select a dressing that is of a suitable size and shape. Ensure the dressing has a 2cm overlap around the wound edge

Dressing Application:

1. Peel the liner from the dressing to expose the adhesive surface
2. Position the dressing centrally over the wound
3. Apply the adhesive side of the dressing onto the skin
4. Gently remove the frame and smooth dressing edges
5. Smooth the entire dressing from the centre toward the edges

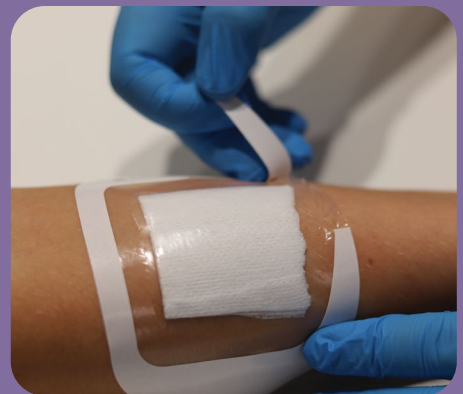

Call: 0800 008 6646

E-mail: woundcare@rociallehealthcare.com

Visit: www.emperorwound.com

 **Rociale**
Healthcare

Emperor is a product brand of Rociale Healthcare